


T-Pro | BEYOND PERFECTION

Material : Gold NiTi Length : 21mm,25mm,31mm Rotary speed : 350 rpm Package : box/blister



Cross Section :

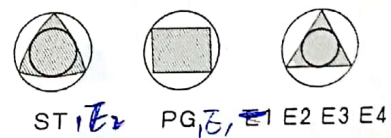

Protocol for use

- Explore the canal using #8/10 or #15 K file, always irrigate the canal when needed.
- Widen the orifice of the canal by ST file.
- Guild path preparation with PG file for curved canal if necessary.
- Remove and sharp with system files E1 & E2 to the working length.
- For severely curved canals, size E3 & E4 are required.

Selling point

- New general files combine Heat-treatment files with Heat-active.
- 200% improved in the durable, enduring cutting force and fracture resistance.
- Files ST/PG/E1/E2 are processed with Heat-treatment while E3/E4 with Pre-bendability.
- Superior quality and best cost-effective for you!

